


EasyFlo™ Uroflow System

EasyFlo Flowmeter
Ver. 2.0 Scale Ver. 1.0
Test ID # 000091

Site #1


Patient: _____
Date: _____ Time: _____

Dr: _____

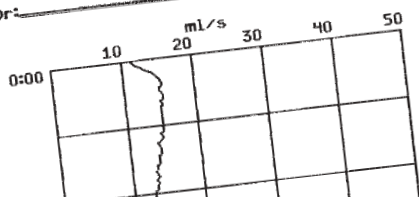

*Accurate screening
for voiding
dysfunction in
any clinical setting*

- ◆ Very simple to operate, with both manual and automatic modes
- ◆ Compact, portable and reliable
- ◆ Full range of clinical measures, including Liverpool nomogram

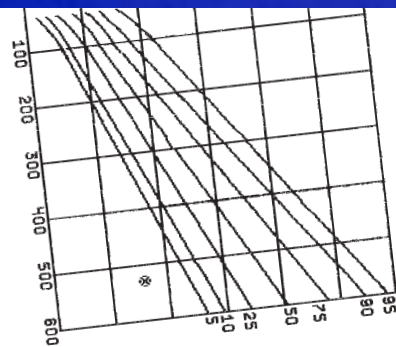

SRS
MEDICAL

Innovative Solutions for Incontinence®

EasyFlo™ Uroflow System

Uroflowmetry remains the mainstay in a urologic workup and is broadly considered the first step in the diagnosis of lower urinary tract symptoms. EasyFlo is designed to meet physician and patient needs for an accurate, easy-to-use uroflow system.

EasyFlo is a compact, portable system that can be used in any clinical setting. Its load-cell technology ensures accurate measurements and, unlike alternative uroflow systems, remains easy to clean. Patient studies print on an easily interpreted report, with all clinical measures automatically charted, including:

- ◆ Urine flow pattern (graphically)
- ◆ Total voided volume
- ◆ Total void time
- ◆ Urine flow time
- ◆ Average urine flow
- ◆ Maximum (peak) flow
- ◆ Onset-to-peak flow time
- ◆ Liverpool Nomogram (Normative comparisons)

EasyFlo automatically starts the test with the onset of urine flow and automatically ends when a period of 30 seconds elapses with no urine flow. There is no need to have personnel present to switch off the device.



Technical Specifications

- ◆ Configuration: Portable
- ◆ Power: Universal input power supply; EN60601-1 marked
- ◆ Selectable Flow Rate Scaling: 1-25 ml/sec or 1-50 ml/sec
- ◆ User's Guide Languages: English, Spanish, French, German, Italian, Dutch, Portuguese
- ◆ Warranty: One year parts and labor


Additional Options

- ◆ Standard Commode Chair
- ◆ Wall Mount/ Tabletop Shelf Assembly
- ◆ Carry Case

Related Products for Continence Care

SRS Medical is a worldwide leader in innovative medical devices for the diagnosis and conservative treatment of incontinence, a condition that affects the quality-of-life of over 30,000,000 adults. In addition to uroflow, SRS Medical offers urodynamics, clinical biofeedback, neurostimulation and the UroScope™ bulking agent injector/female endoscope. SRS Medical patient products include StepFree™ vaginal weights, C3® penile clamp and the In-Flow™ device, an advanced technology for women with atonic bladder (investigational device - not available in U.S.).

- ◆ Compact, portable and reliable
- ◆ Fully automatic or "one-touch" manual operation
- ◆ Automatic report generation with Liverpool nomogram
- ◆ Includes quiet, maintenance-free printer
- ◆ Economical to purchase and to operate

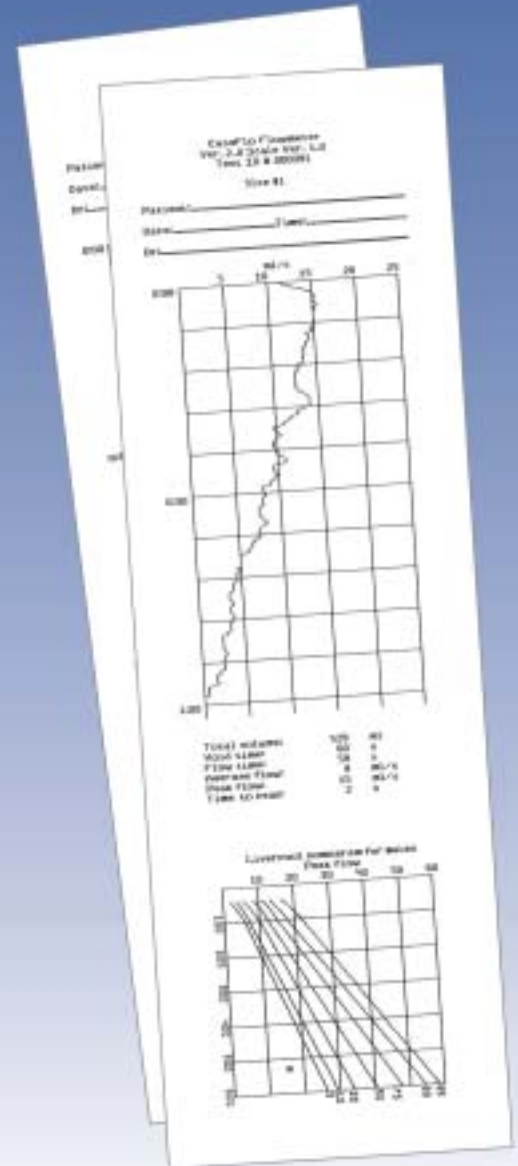

© 2001-2004 SRS Medical Corp. All rights reserved.

Innovative Solutions for Incontinence and C3 are registered trademarks and EasyFlo, UroScope and StepFree are trademarks of SRS Medical Corp. In-Flow is a trademark of Personal Medical Corp.



P/N 21000483 Rev B, 02/04



Innovative Solutions for Incontinence®

800 345 5642 ◆ srsmedical.com